



# FAMED

Żywiec

Member of REINSBERG® GROUP




## Famed VIDAR

SIMPLE OPERATING TABLE  
WITHOUT ELECTRIC CONTROL

# Famed VIDAR

## OPERATING AND TREATMENT TABLE WITHOUT ELECTRIC CONTROL

Functional and affordable operating and treatment table designed to perform uncomplicated procedures for dressing, surgery, orthopaedics, proctology, urology, gynecology, laryngology and others. Simple construction made of high quality carbon steel covered with powder varnish facilitates cleaning and disinfection. The lack of electric control is perfect for the most extreme conditions. Height adjustment is controlled by a hydraulic pump, and other functions are supported by gas springs.

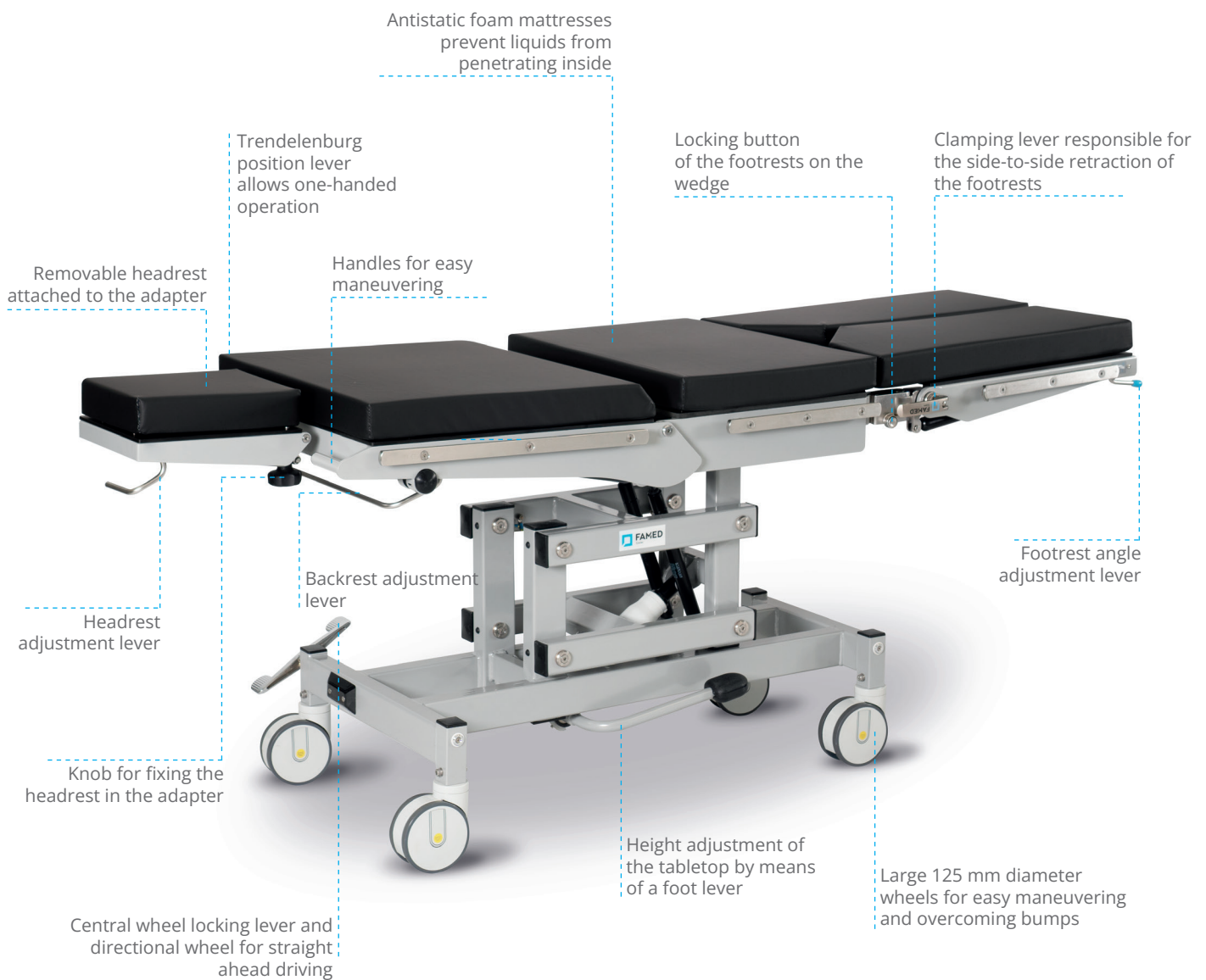

## Key features:

- Pantograph design;
- Antistatic 80 mm thick foam mattresses with waterproof cover to prevent liquid penetration;
- Removable mattresses mounted with hook-and-loop fasteners;
- Table top made of translucent HPL panels;
- Possibility to remove the footrests and headrest and replace them with other accessories;
- Tilt adjustment of footrests, backrest and headrest supported by gas springs;
- A hydraulic pump activated by a foot lever is responsible for adjusting the height of the table top;
- Stable wedge connections with locking mechanism guarantee safe installation of footrests and other accessories;
- Stable base with large wheels makes it easy to move and overcome unevenness;
- The central wheel lock guarantees stability, and the directional wheel makes it easy to go straight;
- Trendelenburg position lever allows one-handed operation;
- The clamping lever allows the footrests to be retracted laterally and placed in the desired position;
- Ergonomic handles located on the side of the patient's head make the table easy to maneuver;
- Full-length side rails allow mounting of various accessories;
- The table design allows patient monitoring and X-rays to be taken during procedures and operations.

## Technical data:

	Famed VIDAR SZ-01
Total length of the table top	2030 mm
Total width of table	650 mm
Width of table top	605 mm
Minimum height of the table top	810 mm
Maximum height of the table top	1110 mm
Backrest adjustment	-40° do +85°
Vertical adjustment of the footrests	-90° do +30°
Horizontal adjustment of the footrests	180°
Headrest adjustment	-45° do +60°
Trendelenburg / anty-Trendelenburg	25° / 15°
Permissible load	150 kg
Weight of table without equipment	125 kg





Famed AS-69 hand surgery tabletop



Arm support Famed FlexArm™ AS-70.3



Stirrups K-TEK



Knee support Famed WS-05.5



Anty-Trendelenburg position 15°



Famed AS-88 knee arthroscopy attachment

Edition 282/2024/03/3 Famed Żywiec reserves the right to modify the product and specification as part of technical progress. All illustrations and photos used in this material are for illustrative purposes only and may not reflect the finished product. The people shown in the photos are not medical professionals. They are models. The device presented in the catalog is intended for use in health care facilities by authorized persons after having read the manual (IFU).

